

V.A.C. Veraflo Cleanse Choice™ Dressing for use with V.A.C. Veraflo™ Therapy



References

1. Kim PJ, Applewhite A, Dardano AN *et al.* Use of a novel foam dressing with Negative Pressure Wound Therapy and Instillation: Recommendations and clinical experience. *Wounds* 2018;30:S1-S17.
2. Teot L, Boissiere F, Fluieraru S. Novel foam dressing using negative pressure wound therapy with instillation to remove thick exudate. *Int Wound J*2017; 14:842-848.

Footnotes

A. Patient data and photos courtesy of Kimberly D. Hall, DNP, RN, GCNS-BC, CWCN-AP, COCN

B. Patient data and photos courtesy of Lindsey Waddell, RN, MSN, WHNP-BC

Note: Specific indications, limitations, contraindications, warnings, precautions and safety information exist for these products and therapies. Please consult a clinician and product instructions for use prior to application. This information is intended for healthcare professionals only.



For more information, visit [Solventum.com](https://www.solventum.com).

V.A.C. Veraflo Cleanse Choice™ Dressing

- When used in conjunction with V.A.C. Veraflo™ Therapy, the V.A.C. Veraflo Cleanse Choice™ Dressing can help facilitate the removal of wound exudate and infectious material such as thick fibrinous exudate and slough.
- Ideal for wound cleansing when surgical debridement must be delayed or is not possible or appropriate.

Why wait weeks when you could see results in days?

V.A.C. Veraflo Cleanse Choice™ Dressing with V.A.C. Veraflo™ Therapy

Can facilitate the following with regards to thick wound exudates such as non-viable tissue:²

Soften → Separate¹ → Solubilize

A consensus panel of surgeon experts provided the following examples of clinical situations in which V.A.C. Veraflo Cleanse Choice™ Dressing with V.A.C. Veraflo™ Therapy could be used in conjunction with appropriate wound care such as debridement and systemic antibiotics: '


- Presence of fibrin, slough, or nonviable tissue.
- Surgical debridement may be delayed or not possible.
- Bed side debridement can not be tolerated.
- Patient chooses not to have surgery.

V.A.C. Veraflo Cleanse Choice™ Dressing recommended wound characteristics:¹


- Majority area of nonviable tissue and have heavy bioburden, and/or are difficult to granulate
- Wounds with exposed bone wounds with treated, underlying osteomyelitis
- Wounds in the presence of orthopedic fixation hardware
- Wounds being prepared for definitive closure or coverage (eg, split-thickness skin graft [STSG], full-thickness skin graft [FTSG], or flap)

Goals for using V.A.C. Veraflo Cleanse Choice™ Dressing are varied and include:

- Cleanse wounds when areas of slough or non-viable tissue remain present on the wound surface
- Remove thick exudates and remove infectious materials
- Promote granulation tissue formation
- Help provide a bridge to a defined endpoint for a clinical plan of care

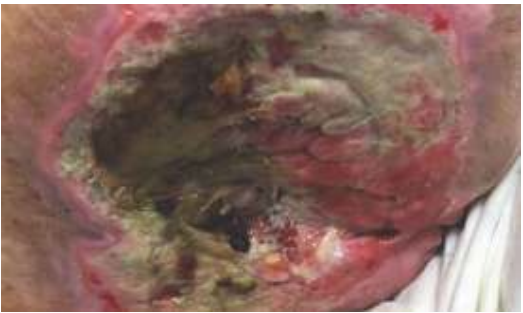


Day 1: The wound had been previously treated with V.A.C.® Therapy, offloading, silver dressings, air mattress use, hydrofiber dressings, alginate dressings, and wound debridement. Bedside sharp debridement was performed but limited by inability to achieve adequate hemostasis.*

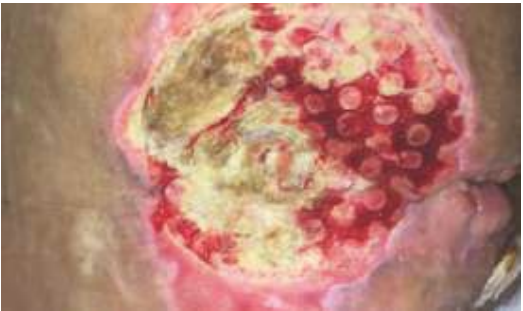

Day 3: V.A.C. Veraflo™ Therapy, using V.A.C. Veraflo Cleanse Choice™ Dressing was initiated. Saline (22 ml) was instilled into the wound followed by a 1 minute dwell time and 30 minutes of negative pressure at -150 mmHg. Due to the difficult wound location, ostomy paste was used to help ensure a complete seal around the wound. After 3 days of therapy, the wound showed improvement.^

V.A.C. Veraflo Cleanse Choice™ dressing with V.A.C. Veraflo™ Therapy proposed mechanism of action

Proposed mechanisms of action of reticulated open-cell foam dressing with through holes (ROCF-CC) combined with negative pressure wound therapy with instillation and dwell time. (A) Noncompressed ROCF-CC dressing with no negative pressure applied; (B) compressed ROCF-CC dressing with negative pressure applied; (C) instillation of topical solution with no negative pressure applied; and (D) wound bed with macrocolumns after dressing removal!

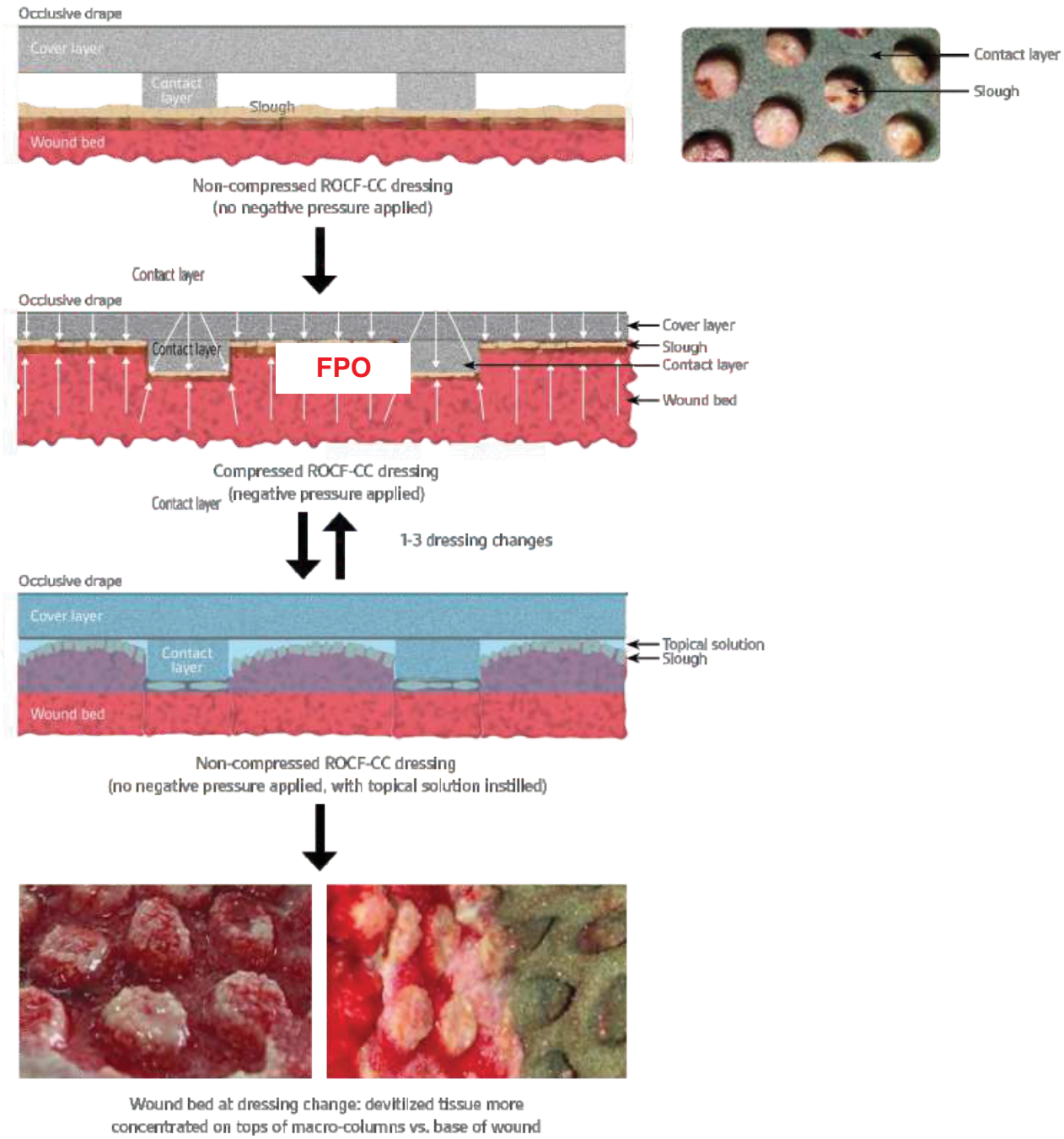

The use of V.A.C. Veraflo Cleanse Choice™ Dressing in a clinical setting

A retrospective data analysis on 21 patients with 21 large complex chronic wounds that contained substantial areas of devitalized tissue and/or yellow fibrinous slough and who were treated in one hospital by several surgeons.

- The V.A.C. Veraflo Cleanse Choice™ Dressing was applied to all wounds using two foam layers: a wound contact layer with 1.0 cm diameter through holes and a cover layer without holes.
 - Dressings were changed every 3 days at the bedside (personal communication). The wounds were washed with saline before a new dressing was applied.
- V.A.C. Veraflo™ Therapy with saline was delivered with the following settings:
 - Soak time: 10 minutes
 - V.A.C.® Therapy phase time: 3.5 hours
 - Target Pressure: -125 mmHg


- Patient Population:

 - 18/21 (85.7%) wounds were pressure ulcers (ishical, sacral and trochanter); 1/21 (4.8%) wounds was a burn wound and 2/21 (9.5%) wounds had necrosis after skin excision.
- 15/21 (71.4%) patients had a confirmed and treated bone infection.
 - 11/21 (52.4%) patients were paraplegic or quadriplegic.

Mean duration of V.A.C. Veraflo™ Therapy with V.A.C. Veraflo Cleanse Choice™ Dressing was 8.7 days with an average of 2.9 dressing changes.

Most of the non-viable tissue was removed at the first dressing change after 3 days of therapy

Within 3 applications (within 9 days) of therapy, wound outcomes observed included:

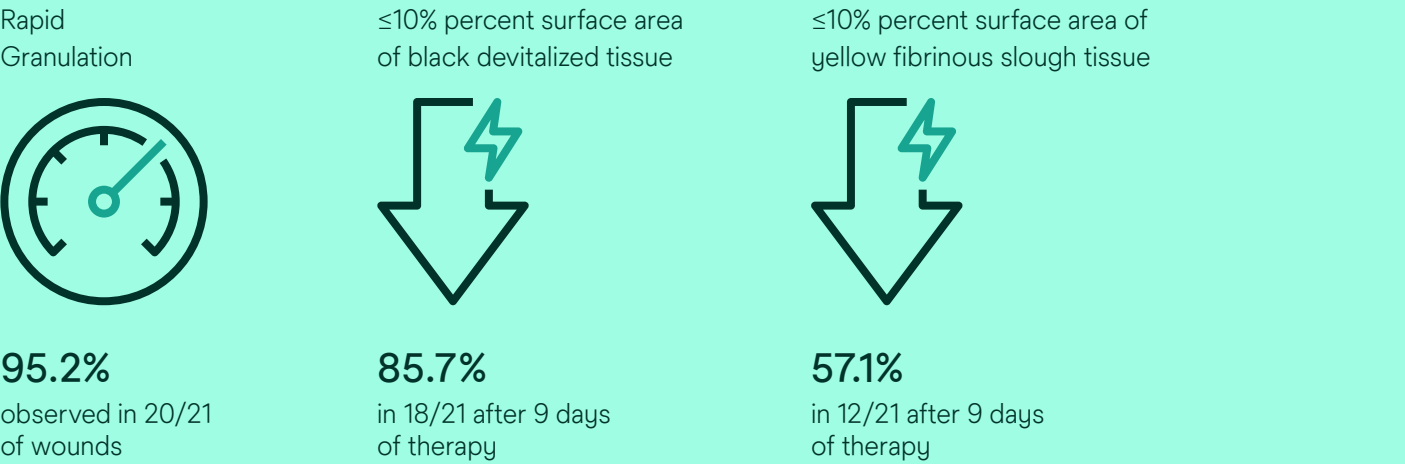


Patient results with V.A.C. Veraflo Cleanse Choice™ Dressing

V.A.C. Veraflo Cleanse Choice™ Dressing with V.A.C Veraflo™ Therapy-stage 4 pressure ulcer of the sacrum^



Day 0: The wound had been previously treated with V.A.C.® Therapy, offloading, silver dressings, air mattress use, hydro-fiber dressings, alginate dressings, and wound debridement. Bedside sharp debridement was performed but limited by inability to achieve adequate hemostasis.

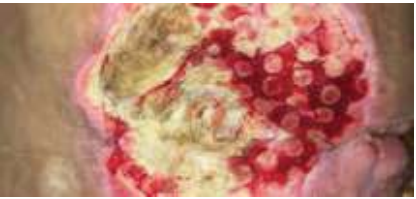

Day 3: The wound showed improvement after V.A.C. Veraflo™ Therapy, using V.A.C. Veraflo Cleanse Choice™ Dressing was initiated. **Saline (22 ml)** was instilled into the wound followed by a **1 minute dwell time** and **30 minutes of negative pressure at -150 mmHg**. Due to the difficult wound location, ostomy paste was used to help ensure a complete seal around the wound.



Day 7: V.A.C. Veraflo™ Therapy was discontinued. The wound underwent sharp debridement to remove the tip of the coccyx and non-viable slough/adipose tissue, followed 2 days later by colostomy surgery. Three days post surgery, V.A.C. Veraflo™ therapy with V.A.C. Veraflo Cleanse Choice™ Dressing was re-started instilling **Saline (22 ml)** followed by a **1 minute dwell time** and **30 minutes of negative pressure at -150 mmHg** (Day 12). A silver alginate dressing was placed over the left buttock partial thickness area.

Patient results with V.A.C. Veraflo Cleanse Choice™ Dressing



Day 0: Venous stasis ulcer measured 22 cm x 16 cm and 85% fibrinous.



Day 10: Wound after negative pressure wound therapy with instillation of **50 ml normal saline** with **3-minute dwell time** followed by **2 hours of -125 mmHg NPWT**.

Patient results with V.A.C. Veraflo Cleanse Choice™ Dressing



Day 0: Flap donor site wound at presentation.



Day 0: wound after bedside debridement.



Day 0: Application of V.A.C. Veraflo Cleanse Choice™ Dressing with V.A.C. Veraflo™ Therapy: instillation of **44 ml of 0.125% hypochlorite solution** with **3-minute dwell time**, followed by **2 hours of -125 mmHg NPWT**.



Day 5: wound appearance.

As with any case study, the results and outcomes should not be interpreted as a guarantee or warranty of similar results. Individual results may vary depending on the patients circumstances and condition