



3M™ RelyX™ Universal
Resin Cement

Goodbye
complexity



The two that can do it all

3M™ RelyX™ Universal Resin Cement works both as a standalone, self-adhesive cement and as an adhesive cement when combined with 3M™ Scotchbond™ Universal Plus Adhesive.

Together they eliminate the hassle and confusion of multiple resin cements, primers and adhesives, **simplifying your direct and indirect restorative workflows**. With fewer products to stock and replenish, you save space, time and money.

So simple yet so powerful



Two components for virtually all indirect and direct bonding indications



Outstanding self-adhesion of RelyX Universal Resin Cement



Scotchbond Universal Plus Adhesive further enhances bond strength to all substrates



Unique 3M automix syringe



Easy excess clean-up



High aesthetics



BPA derivative-free formulations

Game-changing syringe

Proprietary 3M design takes automix cement applications to the next level. Developed with your needs in mind for easy dispensing, cleaning and storage. The small, low-waste 3.4 g syringe provides enough material for approximately 15 to 20 applications.*

Easy excess clean-up

Efficient to the last step. Its unique chemistry ensures excess outflow stays put at the restoration margins, and can be easily removed after tack curing.



Self-adhesive cementation of a 3M™ Lava™ Esthetic Fluorescent Full-Contour Zirconia bridge; Image courtesy of Dr. Jan Gueth, Germany



Less stock, more savings



Easier for operatory staff

95% of dentists agreed the universal character of RelyX Universal Resin Cement helps to **reduce stock and save costs**.

97% of dentists agreed using one universal versus multiple resin cements is **easier for operatory staff**.

* Based upon the amount of cement used in a restoration from 60 µL (e.g. inlay, onlay, veneer) to 90 µL (e.g. crown).

** Per application compared to currently available standard automix systems. Source: M. Cowen, J.M. Powers, Dental Advisor, Biomaterials Research Report, Number of Automix Applications and Mixing Efficiency, No. 134 – June 23, 2020.

*** Source: Field Evaluation EU/USA conducted by 3M.



The strong partner for adhesive cases


3M™ Scotchbond™ Universal Plus Adhesive is the first radiopaque, all-in-one universal adhesive for all direct and indirect bonding indications. And, it is a fully aligned system: the adhesive is cured by the cement, so no light-cure is needed.




Scotchbond Universal Plus Adhesive serves as universal primer for all materials including glass ceramics and further boosts the cement's excellent self-adhesive bond e.g. to zirconia.

Outstanding self-adhesion for maximised efficiency

With its excellent self-adhesive bond strength to dentine, RelyX Universal Resin Cement reliably treats most cases with no need for an additional adhesive or primer. For light- or self-cure.

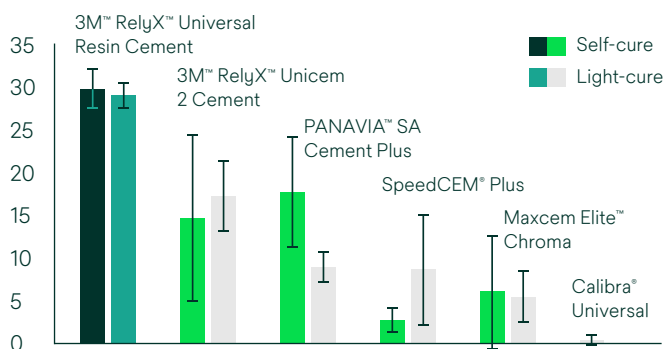

Shear bond strength (MPa) to dentin after artificial aging R. Afutu, M. Abreu, G. Kugel, Tufts University, *J. Dent. Res.* Vol 98A, No 3629, 2019

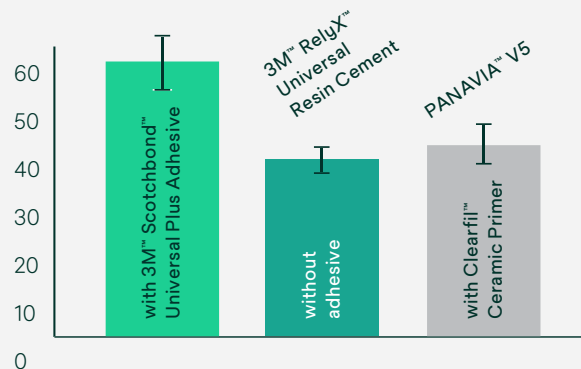

Self-cure shear bond strength (MPa) to zirconia. C. E. Sabrosa et al., *J. Dent. Res.* Vol 99A, No 1838, 2020



Virtually no post-operative sensitivity

99% of dentists reported that RelyX Universal Resin Cement shows virtually no post-operative sensitivity***

High aesthetics with four colour-stable, fluorescent shades



Ordering information

3M™ RelyX™ Universal Resin Cement

Item no.	Product information		
Intro kit			
56968	4 x 3M™ RelyX™ Universal Resin Cement syringes shades TR, A1, A3O, WO (3.4 g each), 60 x 3M™ RelyX™ Universal Micro Mixing Tips, 20 x elongation tips, 1 x vial 3M™ Scotchbond™ Universal Plus Adhesive (5 ml), 1 x 3M™ Scotchbond™ Universal Etchant (3 ml), 25 blue etchant tips, 50 x orange microbrushes, step-by-step guides		
Trial kits			
Each trial kit contains: 1 x 3M™ RelyX™ Universal Resin Cement syringe (3.4 g), 15 x 3M™ RelyX™ Universal Micro Mixing Tips, 5 x elongation tips, 1 x vial 3M™ Scotchbond™ Universal Plus Adhesive (1.5 ml), 1 x 3M™ Scotchbond™ Universal Etchant (3 ml), 25 x blue etchant tips, 50 x orange microbrushes, step-by-step guides			
56969	Shade TR	56970	Shade A1
Refills			
Each refill contains: 1 x 3M™ RelyX™ Universal Resin Cement syringe (3.4 g), 15 x 3M™ RelyX™ Universal Micro Mixing Tips, 5 x elongation tips			
56971	Shade TR	56973	Shade A3O
56972	Shade A1	56974	Shade WO
Value packs			
Each value pack contains: 3 x 3M™ RelyX™ Universal Resin Cement syringes (3.4 g each), 45 x 3M™ RelyX™ Universal Micro Mixing Tips, 15 x elongation tips			
56977	Shade TR	56978	Shade A1
Mixing tips			
56975	30 x 3M™ RelyX™ Universal Micro Mixing Tips		
56976	15 x 3M™ RelyX™ Universal Micro Mixing Tips with elongation tips		



3M™ RelyX™ Universal Resin Cement Intro Kit (56968)



3M™ RelyX™ Universal Resin Cement Trial Kit TR (56969)



3M™ RelyX™ Universal Resin Cement Refill TR (56971)



3M™ Scotchbond™ Universal Plus Adhesive Refill (41294)

3M™ Scotchbond™ Universal Plus Adhesive

Vial			
41293	Intro kit: 1 x 3M™ Scotchbond™ Universal Plus Adhesive vial (5 ml), 50 x disposable brushes, 1 x 3M™ Scotchbond™ Universal Etchant (3 ml), 25 x dispensing tips		
Refill: 3M™ Scotchbond™ Universal Plus Adhesive vial (5 ml)			
41294	1 x 5 ml vial	41295	3 x 5 ml vial
Unit dose			
3M™ Scotchbond™ Universal Plus Adhesive Unit Dose L-Pops, 1 x 3M™ Scotchbond™ Universal Etchant (3 ml), 25 x dispensing tips			
41297	Intro kit: 50 x L-Pop	41299	200 x L-Pop
41298	100 x L-Pop		

3M™ Scotchbond™ Universal Etchant Etching Gel

41263	2 x 3M™ Scotchbond™ Universal Etchant (3 ml each), 50 x dispensing tips
-------	---



Solventum Dental Solutions

KCI Medical Australia Pty Ltd
Level 3, Building A, 1 Rivett North Ryde, NSW 2113
1300 363 454

KCI New Zealand Unlimited
Suite 1701, Level 17, PwC Tower, 15 Customs Street West
Auckland Central, Auckland 1010
0800 808 182

www.solventum.com

Always follow the Instructions for Use (IFU) and refer to IFU for full indications, precautions and warnings. Claims supported by 3M data on file.

© Solventum 2026. Solventum, the S logo, Lava and RelyX are trademarks of Solventum or its affiliates. 3M, the 3M logo and Scotchbond are trademarks of 3M. All other trademarks are owned by their respective owners.