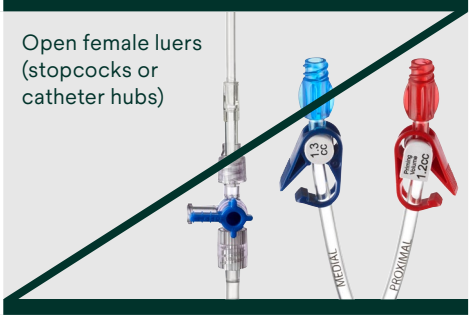
solventum
Curos™

# Cap it with the right colour

| Use this Curos Cap: | Only on: | Do not use on: |
|---|---|---|
| 3M™ Curos™ Disinfecting Cap for Needleless Connectors  | Needleless connectors  | Open female luers (stopcocks or catheter hubs)   |
| 3M™ Curos™ Stopper Disinfecting Cap for Open Female Luers  | Open female luers (stopcocks or catheter hubs)  | Needleless connectors  |
| 3M™ Curos™ Disinfecting Cap for Tego® Haemodialysis Connectors  | Tego® Haemodialysis Connectors  | Open female luers (stopcocks or catheter hubs)   |
| 3M™ Curos Tips™ Disinfecting Cap for Male Luers  | Male luers  | |

solventum

**KCI Medical Limited**
Charnwood Campus
10 Bakewell Road
Loughborough
LE11 5RB

+44 (0)1509 611 611

**3M Healthcare EMEA Export Limited**
Citywest Business Campus
2050 Orchard Avenue
Dublin 24
Ireland

1800 333 377

STERILE R

Active ingredient: 70% Isopropyl Alcohol (IPA). Single use only. This device is not made with natural rubber latex or plasticiser Diethylhexylphthalate (DEHP). Refer to current labeling for indications for use, warnings and product details.

© Solventum 2026. Solventum and the S logo are trademarks of Solventum or its afliates. 3M, the 3M logo and other 3M owned and/or licensed trademarks are trademarks of 3M Company. Other trademarks are the property of their respective owners. OMG2771568.