solventum
Curos™

How to use

3M™ Curos™ Stopper Disinfecting Cap for open female luers

Disinfects in 1 minute. Protects open female luers for up to 7 days.



Peel

- Verify stopcock valve or catheter clamps are in the closed position prior to using the Curos Disinfecting Cap for Open Female Luers
- If open female luer is visibly soiled, clean with 70% Isopropyl Alcohol swab before using cap
- Do not use cap if the seal is broken, torn or punctured
- Detach the cap from the foil
- Once foil is removed, push and twist device clockwise onto the end of the open female luer until tightly connected
- Cap must remain in place for a minimum of 1 minute and may remain on for 7 days if not removed




Push and twist on cathter hub


- The open female luer can be accessed immediately after removing the Curos Stopper disinfecting cap that has been in place for a minimum of 1 minute
- Always place a new Curos Stopper disinfecting cap on open female luer after each use
- Dispose of the Curos Stopper disinfecting cap after every use
- Follow hospital protocol for disinfecting open female luers, if performing sequential accesses

Warning: to avoid potential for injury, use only on open female luers. Ensure the cap is securely tightened.

Caution: potential choking hazard.

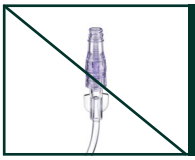

Do not apply to needleless connectors

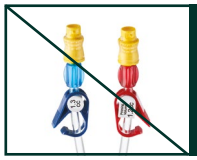

Do not to Tego® Haemodialysis Connectors

Product	Dispenser	Product code	Boxes per case	Units per box	Total caps per case
Curos™ Stopper Disinfecting Caps	Individuals	CSV1-270R	8	270	2,160
	Strip of 10	CSV10-250R	8	25 strips	2,000

solventum

KCI Medical Limited
Charnwood Campus
10 Bakewell Road
Loughborough
LE11 5RB
+44 (0)1509 611 611

3M Healthcare EMEA Export Limited
Citywest Business Campus
2050 Orchard Avenue
Dublin 24
Ireland
1800 33 33 77

STERILE R

Active ingredient: 70% Isopropyl Alcohol (IPA). Single use only. This device is not made with natural rubber latex or plasticiser Diethylhexylphthalate (DEHP). Refer to current labeling for indications for use, warnings and product details.

© Solventum 2026. Solventum and the S logo are trademarks of Solventum or its affiliates. 3M, the 3M logo and other 3M owned and/or licensed trademarks are trademarks of 3M Company. Other trademarks are the property of their respective owners. OMG2771568.