



# A full line of Filtek™ Dental Restorative solutions to fit your clinical needs


With over 2 billion restorations and 50 years of trusted excellence, Filtek™ Dental Restoratives offer reliable solutions with superior esthetics and outstanding durability in virtually every clinical situation.

	Esthetics and versatility		Esthetics and simplicity		Efficiency and simplicity		
Product	<b>3M™ Filtek™ Supreme Ultra Universal Restorative</b> 	<b>3M™ Filtek™ Supreme Flowable Restorative</b> 	<b>3M™ Filtek™ Easy Match Universal Restorative</b> 	<b>Solventum™ Filtek™ Easy Match Flowable Restorative</b> 	<b>3M™ Filtek™ Universal Restorative</b> 	<b>3M™ Filtek™ One Bulk Fill Restorative</b> 	<b>3M™ Filtek™ Bulk Fill Flowable Restorative</b> 
Primary benefits	<ul style="list-style-type: none"> <li>• Exceptional esthetics with natural-looking fluorescence</li> <li>• Outstanding strength for anterior and posterior</li> </ul>	<ul style="list-style-type: none"> <li>• Flows where needed, holds shape without slumping or running</li> <li>• Selection of cannula tip sizes (20G, 27G)</li> <li>• Virtually no bubbles or “run-on”; ergonomic syringe design</li> </ul>	<ul style="list-style-type: none"> <li>• Quick, easy and intuitive shade selection</li> <li>• Naturally-adaptive opacity; no blocker needed</li> </ul>	<ul style="list-style-type: none"> <li>• Quick, easy and intuitive shade selection</li> <li>• Naturally-adaptive opacity; no blocker needed</li> <li>• Virtually no bubbles or “run-on”; ergonomic syringe design</li> <li>• Selection of cannula tip sizes (20G, 27G)</li> </ul>	<ul style="list-style-type: none"> <li>• Streamlined shade offering</li> <li>• Pink Opaquer (PO) to mask dark areas</li> </ul>	<ul style="list-style-type: none"> <li>• One-step bulk placement without sacrificing esthetics and opacity</li> <li>• Increased opacity for improved esthetics</li> </ul>	<ul style="list-style-type: none"> <li>• One-step bulk placement</li> <li>• Virtually no bubbles or “run-on”; ergonomic syringe design</li> </ul>
Maximum incremental placement depth <sup>1</sup>	2mm	2mm	2mm	2mm	2mm	5mm	4mm
Number of shades	36 (Dentin, Body, Enamel, Translucent)	12	3	3	8 (plus Extra White and Pink Opaquer)	5	4
True nanotechnology	Yes	Yes	Yes	Yes	Yes	Yes	No
Handling	Firm	Flows on demand	Firm	Flows on demand	Creamy	Sculptable	Adapts easily
Polish retention & wear resistance	Excellent	Excellent	Excellent	Excellent	Excellent	Excellent	N/A <sup>2</sup>
Safe for composite warming & injection molding <sup>3,4</sup>	Capsules	Syringes and Capsules	Capsules	Syringes	Capsules	Capsules	Syringes and Capsules

For a full list of indications, please refer to the product Instructions for Use.

1. Consult with the instruction for use for the increment depth for each indication. 2. Filtek Bulk Fill Flowable's intended use requires allowing at least 2mm on the occlusal cavosurface for a universal or posterior composite.

3. Warm capsules up to 70°C (158°F) for up to 1 hour. 4. Warm flowable syringes up to 70°C (158°F) for up to 1 hour, up to 25 times.

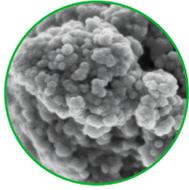

100K magnification of nanocluster. SEM photo courtesy of Dr. J. Perdigão, University of Minnesota.

## Meet the TRUE nanotechnology behind the success of Filtek™ Dental Restoratives

The initial polish of any restorative can be impressive, but many dull over time. Thanks to Solventum's patented true nanofiller technology, Filtek Universal Restorations polish quickly to a high gloss, and polish can even be improved with toothbrush abrasion.\*

\* Solventum internal data

## A family of safe-to-warm composites

From easier handling to reduced extrusion force, composite warming can help your procedure flow more smoothly. Built on a history of clinical success, select Filtek Dental Restoratives have been tested and indicated for composite warming.<sup>1,2</sup>

Now, the portable Solventum™ Filtek™ Composite Warmer makes it quick and easy to bring the benefits of warming to your practice.



1. Warm capsules up to 70°C (158°F) for up to 1 hour.
2. Warm flowable syringes up to 70°C (158°F) for up to 1 hour, up to 25 times.

	Universals	Bulk Fills	Flowables
	3M™ Filtek™ Supreme Ultra Universal Restorative 3M™ Filtek™ Easy Match Universal Restorative 3M™ Filtek™ Universal Restorative	3M™ Filtek™ One Bulk Fill Restorative	3M™ Filtek™ Supreme Flowable Restorative 3M™ Filtek™ Easy Match Flowable Restorative 3M™ Filtek™ Bulk Fill Flowable Restorative
 Class I	✓	✓ (up to 4 mm)	✓ (Liner/Base)
 Class II	✓	✓ (up to 5 mm)	✓ (Liner/Base)
 Class III	✓		✓
 Class IV	✓		
 Class V	✓		✓
Indicated for warming	✓	✓	✓



1-800-634-2249  
[Solventum.com](http://Solventum.com)

You might notice some of our products still carry the 3M logo as we continue our rebranding process. But our top priority is, and always will be, providing you with the exceptional quality and solutions you count on. We never stop solving for you.

© Solventum 2025. Solventum, the S logo and Filtek are trademarks of Solventum or its affiliates. 3M is a trademark of 3M Company.