



Protect surgical incisions, minimize post-operative risk

Prevena™ Therapy is intended to manage closed surgical incisions and surrounding tissue in patients at risk for developing post-operative complications, including surgical site infection.



Prevena™ 125 Therapy Unit



Prevena Plus™ 125 Therapy Unit



Customizable™ Dressing

Peel & Place™ Dressing - 35cm

Peel & Place™ Dressing - 20cm

Peel & Place™ Dressing - 13cm

Versatility meets security

Prevena™ Incision Management System: Indications for Use

Prevena™ Incision Management System is intended to manage closed surgical incisions and surrounding intact skin in patients at risk for developing postoperative complications, such as infection, by maintaining a closed environment via the application of a negative pressure wound therapy system to the incision. The Prevena™ Incision Dressing skin interface layer with silver reduces microbial colonisation in the fabric.

- Disposable, easy to use negative pressure therapy system
- Visual/audible alerts for patient safety
- Replaceable canister
- Patient can shower with Prevena™ Dressing in place

Prevena™ Incision Dressing skin interface layer:

- Actively wicks fluid away from the incision
- Contains 0.019% ionic silver, which reduces bacterial colonization in the fabric

Ordering information

Item #	Product name	Qty
PRE1001	Prevena™ Peel & Place™ System Kit – 20cm	1
PRE1055	Prevena™ Peel & Place™ Dressing – 20cm	5
PRE1101	Prevena™ Peel & Place™ System Kit – 13cm	1
PRE1155	Prevena™ Peel & Place™ Dressing – 13cm	5
PRE3201	Prevena Plus™ Peel & Place™ 35cm System Kit	1
PRE3255	Prevena Peel & Place™ 35cm Dressings	5

The Prevena™ Incision Management Therapy System clinical outcomes:

More than 70+ published clinical studies in a variety of surgical incisions and patient populations have demonstrated that the Prevena™ Incision Management System is associated with positive clinical outcomes compared to standard post-operative dressings.



Surgical Site Infection (SSI) & Seroma

The Prevena™ Therapy System manages surgical incisions by:




Helping to hold incision edges together



Removing fluids and infectious materials



Acting as a barrier to external contamination

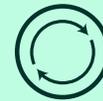

Delivering continuous negative pressure at -125mmHg up to 7 days



Reducing edema

Prevena™ Incision Management System Bibliography available upon request.

Item #	Product name	Qty
PRE4001UK	Prevena Plus™ Customizable™ System Kit	1
PRE4055	Prevena Plus™ Customizable™ Dressing	5
PRE4000	Prevena Plus™ 125 Therapy Unit	1
PRE1095	Prevena™ 45ml Canister	5
PRE4095	Prevena PLUS™ 150ml Canister	5
PRE9090	Prevena™ V.A.C.® Connector	10



For more information, contact your local representative

NOTE: Specific indications, contraindications, warnings, precautions and safety information exist for these products and therapies. Please consult a clinician and product instructions for use prior to application. This material is intended for healthcare professionals.