



Craft smiles with precision and performance

CLARITY™ ADVANCED
CERAMIC BRACKETS





Performance and appearance that exceed expectations

Through the combination of progressive technology and state-of-the-art design, Clarity Advanced Ceramic Brackets deliver brilliant esthetics, small design, predictable debonding and enhanced patient comfort. And, they're available with unique Solventum[™] APC[™] Flash-Free Adhesive, to make the bonding experience faster and more efficient than ever before.



Distinct design

Let your patient's smile shine through in every way with Clarity Advanced Ceramic Brackets. These ultra-small translucent brackets are designed to blend with one's natural tooth color, while resisting staining and discoloration throughout treatment.


Small bracket — strong performance

Showcasing a fine-grained ceramic material and created through an injection-molding process, Clarity Advanced Ceramic Brackets offer trusted strength in a small size and one of the lowest profiles in the industry. Additionally, Clarity Advanced Ceramic Brackets offer ample under tie-wing space, allowing for flexible treatment and ligation options.




Small grain size provides trusted strength in a smaller bracket

*Equal magnification.
Source: Solventum lab.*

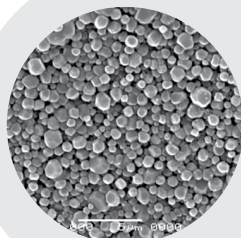


Fine-grained polycrystalline ceramic material—Clarity Advanced Ceramic Brackets

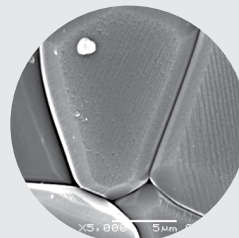

Coarse-grained polycrystalline ceramic material



Clarity Advanced Ceramic Brackets

A revolution in ceramic bracket technology

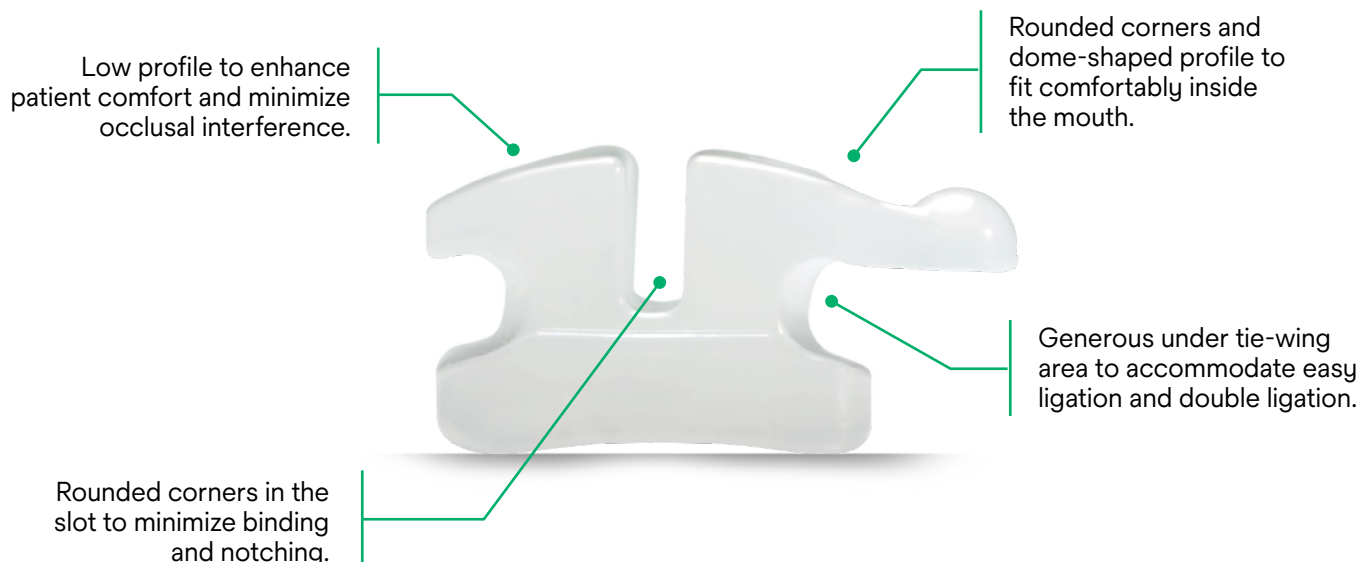

Treat esthetically without compromise

Available in the MBT™ Appliance System,
Roth* and High Torque prescriptions.

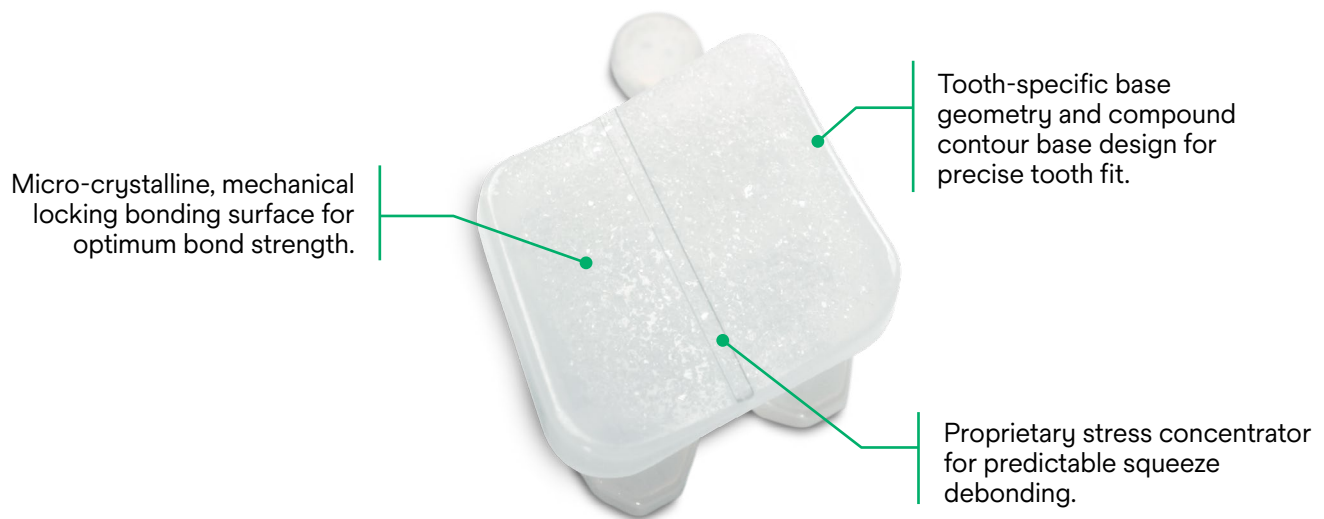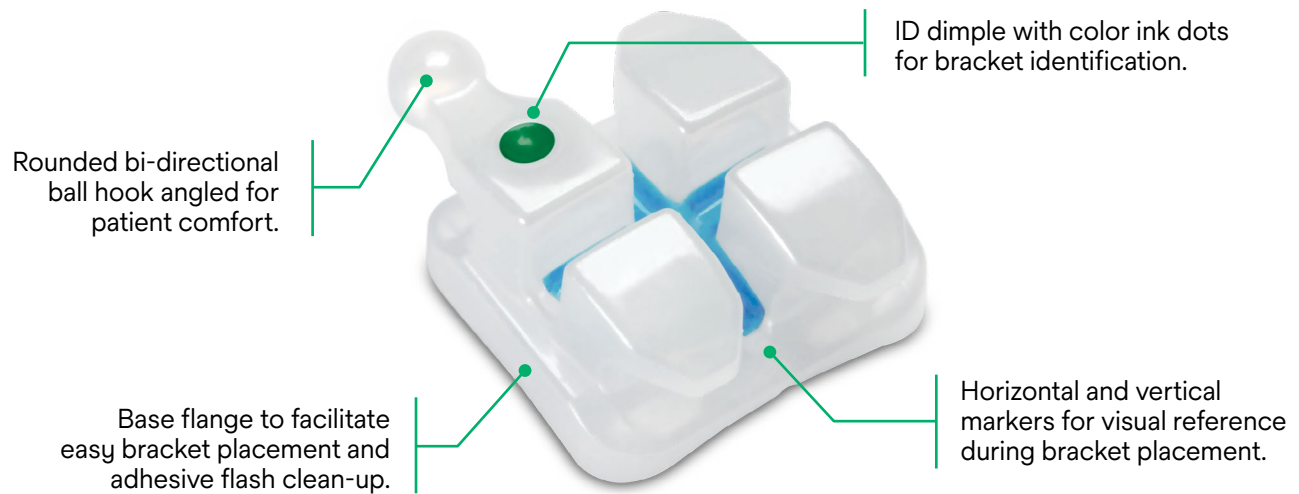


Clarity Advanced
Ceramic Brackets



Forsus™ Class II
Correctors



Tooth-colored
ligature ties

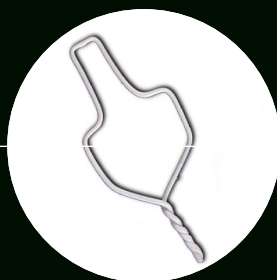



Average bracket
bond failure rate of
less than

2%*



Available precoated
with APC Flash-Free
Adhesive for superior
bonding confidence

Simple, efficient bonding process

Powered by technology, the APC Flash-Free Adhesive Coated Appliance System is the most efficient bonding system in orthodontics. Our innovative design creates a unique adhesive bonding substrate on the base of the bracket that saves you time and positively impacts patient satisfaction.

Saves time and cost

- ✓ Faster and more efficient bonding than non-precoated brackets
- ✓ Eliminates the flash removal step from bonding
- ✓ Simple, easy clean-up at debonding

Superior bonding efficiency

- ✓ Bond time is reduced up to 40% per bracket
- ✓ No flash removal reduces accidental bracket movement after positioning, reducing the need to reposition
- ✓ Protects tooth enamel under adhesive

*within 90 days of bonding

Designed for happy patients and efficient practices

Predictable, consistent debonding with proprietary stress concentrator

Reduce the uncertainty of the ceramic bracket debonding process with our unique stress concentrator, a hallmark of the Clarity brand. This innovative design feature allows for mesial-distal debonding with a simple squeeze, so you can debond with confidence and help ensure the patients' debonding visit will be a positive part of the treatment experience.



Blend in *or* be bold

Patients can choose from more than twenty clear or colorful AlastiK™ Easy-to-Tie Ligatures for any occasion. Small, low-profile bracket with smooth bracket surface and rounded corners enhance patient comfort.

Practice with clarity

Our full portfolio of naturally esthetic treatment options are informed by science and fueled by creativity. Innovation drives our leading-edge technology and delivers the choice and control you seek.

Our commitment is to collaborate with orthodontists to provide custom, flexible esthetic treatment solutions that showcase your expertise and ensure precise outcomes. Together, we can change lives by creating healthy, beautiful smiles.

Clarity™ custom solutions



Clarity™ Ultra Self-Ligating Brackets

Treat your patients to better experiences with innovative, fully-esthetic ceramic brackets designed for patient comfort and reliable results throughout treatment.



Clarity™ Aligners

Embrace the power of customization with a system of removable, nearly invisible, clear aligners that deliver control, flexibility and choice.



Clarity™ Advanced Ceramic Brackets

Enjoy benefits both you and your patients have been looking for — brilliant esthetics, small design, predictable debonding and enhanced comfort.



Solventum Dental Solutions
2510 Conway Avenue
St. Paul, MN 55144-1000 USA
Phone 1-800-423-4588
Web [solventum.com](https://www.solventum.com)

Solventum Canada
300 Tartan Dr.
London, ON N5V 4M9 Canada
Phone 1-800-443-1661

© Solventum 2025. All rights reserved.
APC, Clarity, Forsus and MBT are
trademarks of Solventum. Used under
license in Canada.